

Marlen UltraMax™

One-Piece Disposable Ostomy System



MAXimum Flexibility

MAXimum Comfort



MAXimum Security



*For MAXimum Choice And
MAXimum Flexibility In A One-Piece System* **Marlen UltraMax™**




**UltraMax™ Drainable Pouch For
Ileostomies and Colostomies**

- Extra-large drainage spout for easy emptying
- Marlen Kwick-Klose™ Fastener**
- A secure, dependable closure for drainable pouches
 - Hook & loop closure eliminates the need for a separate clamp
 - Bottom of pouch opens easily for draining


**UltraMax™
Features**

- Advantages of a one-piece system
- Ease of a disposable system
- Odor-proof pouch material
- Quiet, built-in comfort cover facing the body
- Leak-proof for added protection and security

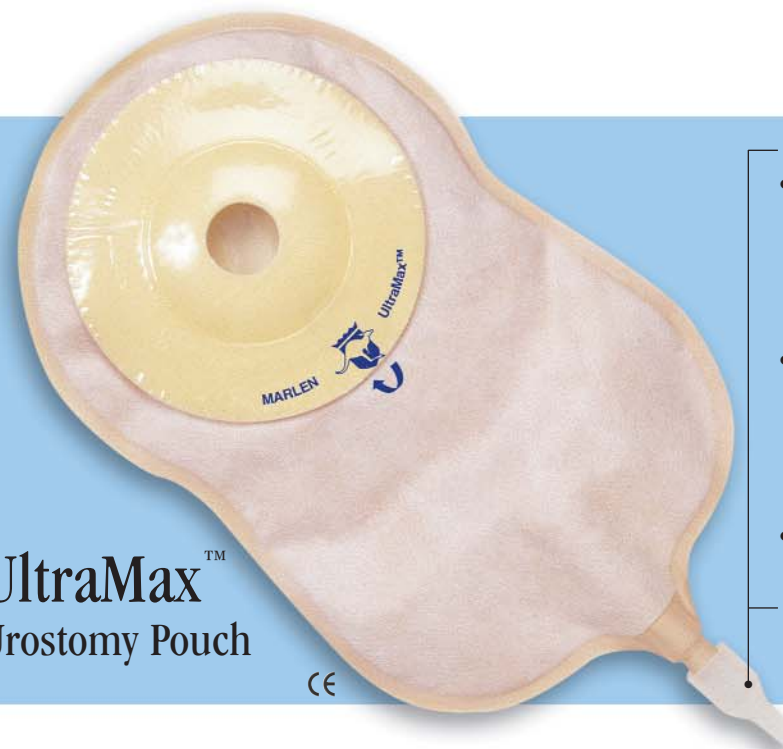

**UltraMax™
Urostomy Pouch**

- “Free-Flow” stem allows unobstructed urine flow into leg or night bags, even if the pouch twists
- “Push-Pull” E-Z Drain Valve is permanently sealed to the pouch. No wet fingers when emptying pouch.
- Non-Return Anti-Reflux Valve

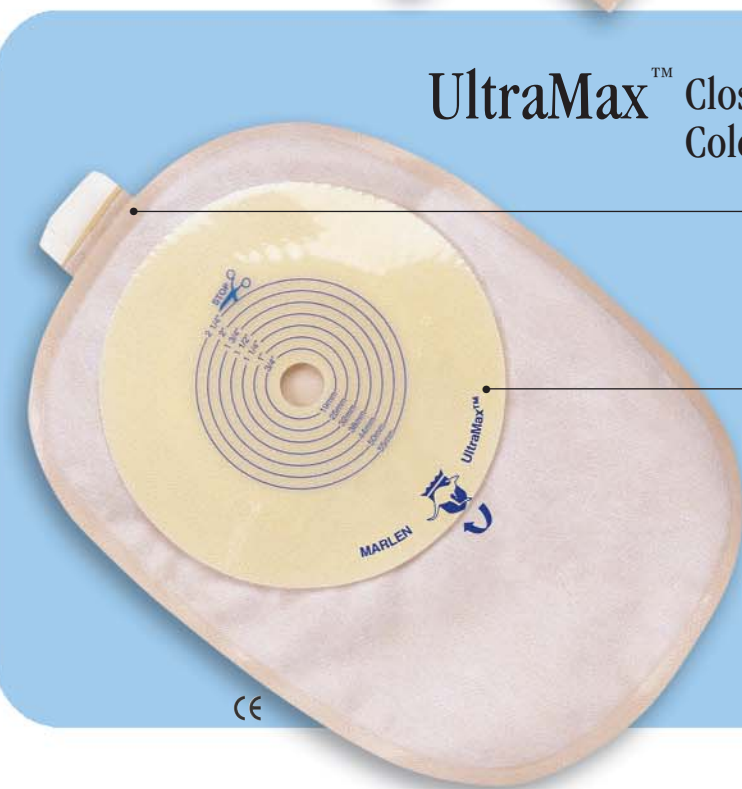

**UltraMax™ Closed End Pouch For
Colostomies**

- “Flo-Thru™” Charcoal Filter System continuously filters odors from both sides of the filter thanks to its patented central suspension design.
- Large adhesive area for secure fit



**AquaTack™
Hydrocolloid Barrier**

- Covers entire flange for maximum skin protection
- Absorbs skin moisture
- Adheres to moist skin
- Seals securely and removes easily from sensitive skin (controlled adhesion)

**Cut-To-Fit
Shallow
Convexity**



Cut-To-Fit design for precise fit and comfort, even for oval and irregular shaped stomas
Starter hole is ½" (12mm). Cut to fit stoma sizes up to 1½" (38mm).

**Cut-To-Fit
Flat**



Cut-To-Fit design for precise fit and comfort, even for oval and irregular shaped stomas
Starter hole is ½" (12mm). Cut to fit stoma sizes up to 2¼" (55mm).

**Pre-Cut
Shallow
Convexity**



Pre-Cut stoma openings available from ½" (12mm) to 1¾" (44mm) in ¼" (3mm) increments.

**Pre-Cut
Flat**




Pre-Cut stoma openings available from ½" (12mm) to 2¼" (55mm) in ¼" (3mm) increments.

MAXimum Flexibility


UltraMax™ Specifications

UltraMax™ Closed End

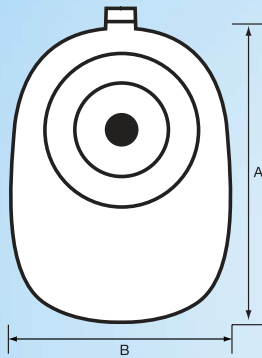


Length A	Width B	Capacity (approx.)
7-1/2" (19 cm)	5-1/2" (14 cm)	9 oz. (270 cc)

UltraMax™ Drainable

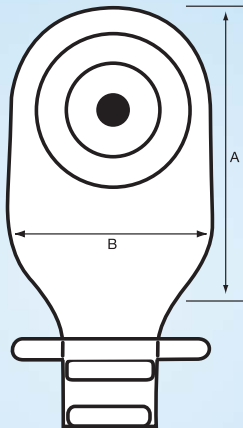

Length A	Width B	Capacity (approx.)
9" (22.9 cm)	5-3/4" (14.6 cm)	16 oz. (480 cc)

UltraMax™ Urostomy

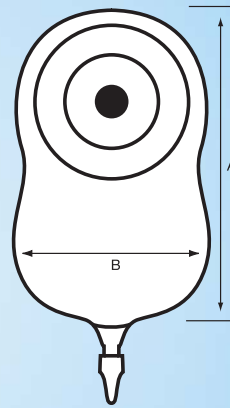

	Length A	Width B	Capacity (approx.)
Regular	8-3/4" (22.2 cm)	5" (12.7 cm)	14 oz. (420 cc)

MAXimum Comfort

MAXimum Security



Marlen Manufacturing & Development Company
 5150 Richmond Road • Bedford, OH 44146 USA
 Phone: 216.292.7060 • Fax: 216.292.9196
 E-mail: info@marlenmfg.com • www.marlenmfg.com